

## HEMOSTAT DRESSING/INTRANASAL SPLINT




- Biocompatible-made of all natural Chitosan-based polymers
- Minimizes bleeding and edema after surgery
- Easily conforms to nasal anatomy
- Offers easier placement with a superior field of view
- No special storage conditions required

## PRODUCT DESCRIPTION:

**PosiSep®** is a unique carboxymethylated chitosan-based hemostatic intranasal splint that has been custom formulated for OR and Office Based FESS procedures with minimal bleeding and where rapid clearance of the product is desirable. After hydration, **PosiSep®** will immediately expand, turn into a gel, then begin clearing the sinus cavity. Subsequent irrigation will cause the product to dissolve, thus optimizing patient tolerance and comfort.

- Dissolves naturally by way of gentle irrigation and aspiration
- Easy to apply-no mixing or additional prep time necessary

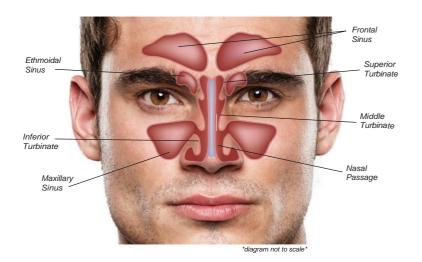


**PosiSep®** Hemostat Dressing/Intranasal Splint is indicated for use in patients undergoing nasal/sinus surgery as a space occupying splint and hemostat to:

- Separate tissue or structures compromised by surgical trauma;
- Separate and prevent adhesions between mucosal surfaces during mesothelial cell regeneration in the nasal cavity;
- Help control minimal bleeding following surgery or trauma;
  Help control minimal bleeding following surgery or nasal trauma by tamponade effect, blood absorption and platelet aggregation; and
- Act as an adjunct to aid in the natural healing process.

PosiSep is indicated for use as a nasal hemostat to treat epistaxis. PosiSep is intended for use under the direction of a licensed healthcare provider.



## Before Hydration After Hydration

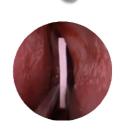



When PosiSep® is hydrated, its smooth gel action separates sinus mucosal tissue.

## What is Chitosan?

shells of crustaceans.

Chitosan (ky-toe-san) comes from a substance known as chitin, an amino-polysaccharide removed from the powdered



Norsk Medikal as Postboks 456 5501 HAUGESUND Tlf. 5270 4140 Fax 5270 4144 mail@norsk-medikal.no www.norsk-medikal.no

ORDERING INFORMATION PosiSep® is available in two sizes. See ordering information below.

Part #: Description: How Supplied:

9215654 PosiSep® Hemostat Dressing/Intranasal Splint 1.5 cm x 5.1 cm/ 8pk
 9215644 PosiSep® Hemostat Dressing/Intranasal Splint 1.3 cm x 3.8 cm/ 8pk

PosiSep® is packaged eight (8) sterile, single use units per carton. One (1) sterile unit contains one (1) Hemostat Dressing/Intranasal Splint. PosiSep® is terminally sterilized by gamma irradiation.

